



PRODUCT CATALOG

MedSafety Solutions is committed to providing you with the pharmacy supplies you need, when you need them – and we've put the people, products and processes in place to make that happen.

As you review and select products from this catalog, you can rest assured that our friendly, knowledgeable team will play our part in keeping your pharmacy running smoothly. If you have any questions about selecting, ordering or using our products, we stand ready to get you the information you need and the service you expect.

Your Job is to Protect and Save Lives; Our Job is to Help You Do It.



KwikDose™

Oral Dispensers



KwikDose oral dispensers accurately deliver small volumes of oral liquid medication. Enlarged numbers and easy-to-read marking lines allow for quick and accurate volume identification. The unique tip prevents wrong-route administration.

FEATURES

- High-quality plunger easily draws and dispenses oral medication
- Large, legible lines provide clear dosage information
- Six convenient sizes: 0.5, 1, 3, 5, 10, and 20 mL
- Available in both clear or amber
- Oral tip designed to prevent wrong-route medication delivery
- Unique, self-righting tip cap

INDICATIONS FOR USE

- Oral dispenser is intended to administer a measured amount of oral liquid medication.
- Oral tip cap is for use with the oral dispenser to reduce the risk of touch contamination and medication leakage.

SPECIFICATIONS

- Non-sterile
- ISO 80369-1 compliant
- Meets ISO 7886-1 accuracy requirements
- Not made with natural rubber latex
- Oral Dispenser - Made in Denmark
- Oral Tip Cap - Made in the USA
- Packaged in the USA

PRECAUTIONS


- Oral dispensers are for oral use only.
- Verify compatibility between drugs and drug contact materials of this product with the drug manufacturer's package insert and/or material safety documentation.
- Do not reuse oral dispensers or tip caps. Product reuse can be unsafe because of the risk of infection, cross contamination, endotoxin reaction, patient injury due to material failure and/or chemical reaction due to cleaning or sterilization agents.
- Dispose of properly according to applicable disposal rules and regulations.
- Rapid dispensing of liquid into mouth may cause choking.
- Oral tip caps are a choking hazard and should be kept out of the reach of children.

ORDERING INFORMATION

mL		MEDSAFETY	ABC/CENCORA	CARDINAL	MCKESSON	MORRIS & DICKSON	HEALTH CARE LOGISTICS
CLEAR ORAL DISPENSER							
0.5	100	7100	10296654	5973797	–	–	21184C
0.5	500	7500	–	–	–	–	21184C-31
1	100	7101	10297801	5973805	–	–	21185C
1	500	7501	–	–	–	–	21185C-31
1	3000	74201	–	–	–	–	–
3	100	7103	10297802	5973813	–	–	21186C
3	500	7503	10302433	–	–	–	21186C-31
3	2500	74203	–	–	–	–	–
5	100	7105	10297660	5973821	3024882	–	21187C
5	500	7505	–	–	–	–	21187C-31
5	1500	74205	–	–	–	–	–
10	100	7110	10297601	5973839	–	–	21188C
10	500	7510	–	–	–	–	21188C-31
10	900	74210	–	–	–	–	21189C
20	100	7120	10297459	5973847	–	–	21189C
20	500	–	–	–	–	–	21189C-31
AMBER ORAL DISPENSER							
0.5	100	8100	10297457	5973938	3033099	–	21184A
0.5	500	8500	10302150	–	–	–	21184A-31
1	100	8101	10297595	5973946	3033107	–	21185A
1	500	8501	10302065	–	3027794	–	21185A-31
1	3000	84201	–	–	–	–	–
3	100	8103	10297630	5973953	3033115	–	21186A
3	500	8503	10302035	–	3024098	–	21186A-31
3	2500	84203	–	–	–	–	–
5	100	8105	10297573	5973961	3027810	–	21187A
5	500	8505	10302068	–	3024114	–	21187A-31
5	1500	84205	–	–	–	–	–
10	100	8110	10297572	5973979	3027828	–	21188A
10	500	8510	10302066	–	3027802	–	21188A-31
10	900	84210	–	–	–	–	–
20	100	8120	10297579	5973987	3027836	–	21189A
20	500	–	–	–	–	–	21189A-31
ORAL TIP CAPS							
–	100	50300	–	–	–	–	–

Contact your distributor for additional ordering information.

KwikVial® Oral Liquid Containers

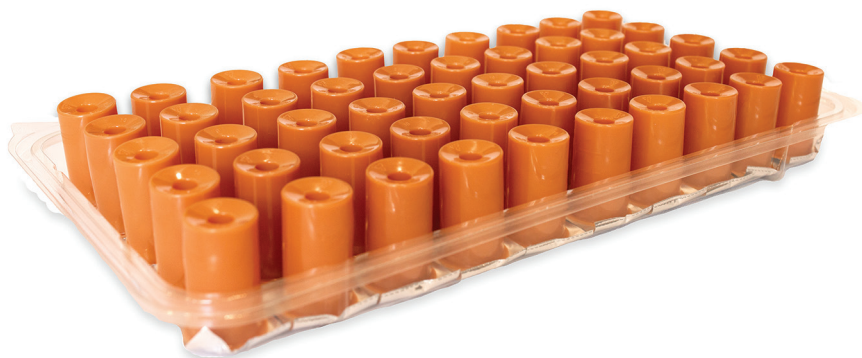

KwikVial Containers are used for unit dose packaging, delivery and administration of oral liquid medications. Versatile and easy-to-use, KwikVial Containers can accommodate a variety of dosage volumes, are sealed for tamper prevention and detection, and are designed for efficient packaging and administration.

FEATURES

- **EFFICIENT** – The vial's bottom-located opening allows the vial to be filled easily by hand or with an automated pump. The KwikVial Filling Extension is designed to connect with any male luer fitting to enable rapid, secure fluid transfers into the vial.
- **SECURE** – The factory-sealed foil lid cannot be opened and re-sealed without detection, minimizing the risk of misuse. The press-in plugs add another layer of protection, creating a tight seal at the container bottom. The plugs can be inserted quickly and efficiently with the fifty-plug capacity KwikVial Plug Dispenser.
- **VERSATILE** – Available in 15 mL and 30 mL versions; both are marked with 5 mL reference graduations, aiding in accurate dosing in 5 mL increments up to 30 mL.
- **SPACE-SAVING** – Designed for efficient storage, two hundred 15 mL container use less than one cubic foot of storage space

INDICATIONS FOR USE

The KwikVial is an amber plastic vial with heat-sealed foil lids, used for unit dose packaging, delivery and administration of oral liquid medications.

SPECIFICATIONS

- Available in two sizes:
 - » 15 mL
 - » 30 mL
- 5 vials per strip connected with easily separable perforations
- 300 vials per case (6-50 vial trays)
- Sealing plugs included
- Translucent amber vial conforms to USP light transmission and tightness requirements.
- Non-sterile
- Made in USA

CONTRAINDICATIONS/PRECAUTIONS

- Verify compatibility between drugs and KwikVial materials with the drug manufacturer's package insert and/or material safety documentation.
- Do not reuse.

KwikVial® Accessories

Make easy-to-use KwikVial Oral Liquid Containers even easier to use and more efficient with these helpful accessories:



PLUG DISPENSER

Enables plugs to be inserted quickly and efficiently, creating a tight seal at the container bottom. The KwikVial Plug Dispenser has a fifty-plug capacity.

- Sold individually
- Non-sterile
- Made in USA



FILLING EXTENSION

Easily attaches to any male luer fitting like pump tube sets, repeating syringes or large fill syringes. The Filling Extension tip extends into filling hole of KwikVial Container, guiding the fluid path and protecting against spillage and waste.

- Clear PVC
- Bulk packaged
- 50 per bag
- Non-sterile
- Made in USA

ORDERING INFORMATION

PRODUCT	MEDSAFETY	ABC/CENCORA	CARDINAL	MCKESSON	MORRIS & DICKSON	HEALTH CARE LOGISTICS
KwikVial Container, 15 mL	30015	10178188	1848548	3656303	894329	7733
KwikVial Container, 30 mL	30030	10178186	1848555	3656329	894337	7735
KwikVial Plug Dispenser	30009	10178187	–	–	–	7339
KwikVial Filling Extension	30008	10178185	1607167	3656337	374207	20226

TwoFer® Needle



TwoFer Huber-Tipped Needles enable vented and non-vented vial additions and withdrawals, effectively performing the functions of two devices. The TwoFer Needle facilitates the movement of liquid drugs from one container to another without the time and expense of changing needles. Color-coded and available in 1.5" and 1.0" lengths, the TwoFer Needle is appropriate for a variety of receptacle sizes and applications.

FEATURES

- **EFFICIENT** – Whether you're using fluid transfer pumps or repeating syringes, the dual-purpose TwoFer Needle enables efficient reconstitution and transfer without changing needles, saving procedural steps and valuable time.
- **EASILY IDENTIFIABLE** – TwoFer hubs are color-coded for quick, easy identification.
- **EFFECTIVE & VERSATILE** – Unique configuration allows for use as a vented or non-vented needle, reducing the handling and extra time required for needle changes.
- **SAFE** – The precision ground stainless Huber point permits smooth and safe vial access.

INDICATIONS FOR USE

Designed as a needle for withdrawing medication from a vial using the metal cannula and for adding to a vial by using the plastic vent hub to equalize pressure.

SPECIFICATIONS

- 2 lengths to choose from:
 - » **LONG** (Green vent/hub)
Approximately 0.5" vented hub plus 1" cannula
 - » **SHORT** (Purple vent/hub)
Approximately 0.5" vented hub plus 0.5" cannula
- 16 gauge
- Silicone coated
- Individually packaged in sterile and non-pyrogenic, undamaged cap and sheath (sterile barrier).
- 100 per case
- Huber point
- Made in the USA


CONTRAINDICATIONS/PRECAUTIONS

- For pharmacy use only, not for use with blood administration.
- Verify compatibility between drugs to be transferred and TwoFer Needle materials with the drug manufacturer's package insert and/or material safety documentation.
- Do not reuse.

ORDERING INFORMATION

PRODUCT	MEDSAFETY	ABC/CENCORA	CARDINAL	MCKESSON	MORRIS & DICKSON	HEALTH CARE LOGISTICS
TwoFer Needle Long/Green	19003	10021369	1503408	3742145	894006	20225
TwoFer Needle Short/Purple	16004	10021390	1887439	3778636	428136	20224

Tamper-Resistant Add Port Cap

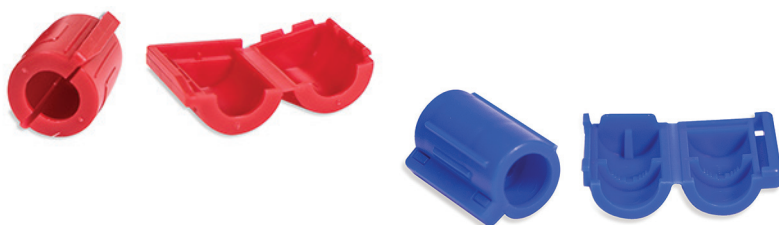

Tamper-Resistant Add Port Caps are easy-to-use, hinged plastic caps that impede tampering when storing and transporting IV solutions.

FEATURES

- **EFFICIENT** – Easy to apply and designed for one-handed application.
- **EASILY IDENTIFIABLE** – Available in red or blue for use with Add Ports of differing sizes (see specifications).
- **EFFECTIVE** – Caps are constructed of thick, durable polypropylene or high-density polyethylene that prevent needle penetration. When applied, the caps are difficult to remove, effectively restricting port access.

INDICATIONS FOR USE

For use on IV bag additive port to prevent tampering by restricting port access.

SPECIFICATIONS

- **RED** – For Add Ports with a Tubing OD less than 0.380", a Shoulder OD between 0.390" – 0.550", and a Shoulder Height less than 0.746". Known to fit Baxter® EXACTAMIX® TPN Containers and Douglas Medical Products EcoFLX™ IV Containers.
- **BLUE** – For Add Ports with a Tubing OD less than 0.34", a Shoulder OD between 0.28" – 0.37", and a Shoulder Height less than 0.37". Known to fit Baxter Viaflex® and similar bags.
- 100 Per Bag
- Non-sterile
- Made in the USA

CONTRAINDICATIONS/PRECAUTIONS

- If the cap cannot be secured (snap closed) or if the additive port diameter is smaller than the cap opening, the Tamper-Resistant Add Port Cap should not be utilized.
- Do not reuse

ORDERING INFORMATION

PRODUCT	MEDSAFETY	ABC/CENCORA	CARDINAL	MCKESSON	MORRIS & DICKSON	HEALTH CARE LOGISTICS
TRAP Cap, Red	MSS4858	10287759	5559091	3962768	489781	7979-03 (25 pc)
						7979-31 (25 x 25 pc)
						7979-BULK (10 x 100 pc)
TRAP Cap, Blue	12020	10287923	5708425	2309144	050021	4073 (25 pc)
						4073-31 (3 x 300 pc)

Female Luer Lock Cap



The Female Luer Lock Cap is designed to seal male Luer fittings on devices such as manifolds, stopcocks or IV sets, reducing the risk of touch contamination and medication leakage.

FEATURES

- **EFFECTIVE** – Reduces the risk of touch contamination and medication leakage.
- **SAFE** – Installs quickly and efficiently; supports aseptic methodology.
- **EFFICIENT** – Simple, fast installation supports a variety of capping requirements.

INDICATIONS FOR USE

For use as a cap for male Luer fittings on medical devices such as manifolds, stopcocks or IV sets.

SPECIFICATIONS

- Individually packaged to be sterile, non-pyrogenic and disposable.
- 100 per case
- Made in the USA

CONTRAINDICATIONS/PRECAUTIONS

- Not intended for blood administration sets.
- The Luer Lock Cap terminates the fluid path.
- Caps are a choking hazard, keep out of reach of children.
- Do not reuse

ORDERING INFORMATION

PRODUCT	MEDSAFETY	ABC/CENCORA	CARDINAL	MCKESSON	MORRIS & DICKSON	HEALTH CARE LOGISTICS
Female Luer Lock Cap, Red	75309	10253827	5708433	2310605	050039	7857-31 (100 pc)
						7857-10 (5 pc)



11233 E. Caley Ave, Suite 300
Centennial, CO 80111

720.500.2500

720.500.2571 (fax)

info@medsafetysolutions.com

www.medsafetysolutions.com